

Therm-X = Heat + Cold + Contrast + Compression

Imagine not having to keep refilling, cleaning and drying out your ice based machine. Imagine being able to maintain the specific temperature you want delivered to the patient, for as long as you want. Then consider the new treatments you can apply with heat, compression, and contrast.

Easy fill tank uses distilled water and isopropyl alcohol, does not require constant refilling and cleaning

Convenient carry handle makes the Therm-X easy to transport at only 14 lbs.

Liquid cooled radiator & fan technology allows the Therm-X garments to get colder and the machine to stay quieter

Umbilical quick connect for a sure and secure connection

Intuitive touch screen with pre-programmed or customized protocols



USB charger for patient's phone/tablet during treatment

NO ICE NEEDED

Therm-X is a Class II Medical Device and should only be operated by a licensed medical professional

Therm-X® uses a new thermal-control technology to treat post-traumatic and post-surgical conditions where localized thermal therapy is indicated.

This **innovative technology** allows the user to select Heat, Cold, or Contrast- temperature range of 34°F to 100°F.

Four garment compression settings ensure full contact for **optimal thermal transfer**.

With its **digital control**, the Therm-X circulates a coolant that can maintain a specific temperature for the entire treatment cycle. Compare this to heat packs or ice-based machine that lose temperature once they are applied.

The fluid passes through the garment with a **tailored treatment** area to concentrate the cold or heat where it is needed the most.

Therm-X® Garments save you money and time

Therm-X® garments are **bilateral** so you don't have to buy lefts and rights. They're **one size fits all**, and have high-end Velcro to last the life of the garment. The inner face uses a ripstop nylon with a waterproof backing to prevent buildup of mold and mildew after cleaning.



Shoulder/XL
Shoulder
TX0101/TX0110



Elbow TX0103



Back TX0105



Therm-X AT Unit
TX0002



Knee TX0102



Hip TX0108



Ankle TX0104



Half Leg Boot TX0111

Accessories to extend and enhance the use of your Therm-X® unit and Garments:



Split Umbilical Hose
TX0208




1 Qt Coolant
TX0206



Extender Straps
TX0109



Carry Case
TX0202

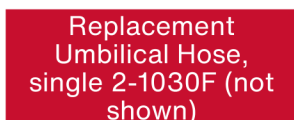

Replacement
Umbilical Hose,
single 2-1030F (not
shown)



Hospital Grade Power
Supply TX0207



Replacement Power
Supply 4-1531